



Corrosion Inhibitor for the Pre-Treatment of Surgical Instruments.



Corrosion inhibitor for spraying onto surgical instruments in operating theatres directly after use

- For the prevention of pitting corrosion caused by residues of saline solution
- Self-acting pre-cleaning already in the OR
- For surgical instruments, incl. MIS instruments
- neodisher

 PreStop is not a hazardous product according to the CLP-Regulation (EC)
 No 1272/2008.




Safe corrosion prevention:

 neodisher
 PreStop prevents pitting corrosion on stainless steel instruments caused by adhesive residues of blood or physiological saline solution.



Application recommendation:

- **neodisher PreStop** is sprayed undiluted directly after surgery onto the surgical instruments in the disposal container ensuring a complete wetting. Keep disposal container closed.
- Instruments treated with neodisher
 PreStop can be stored for up to 72 hours.
- Afterwards rinse carefully and reprocess as usual.

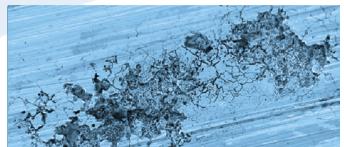

Damage caused by pitting corrosion on stainless steel instruments is irreparable.

Our Competence – Your Benefit:

- reliably prevents pitting corrosion
- keeps the instruments wet
- contains a self-acting pre-cleaning already in the OR
- inhibits the growth of micro-organisms on instruments during contact time (bacteriostatic)
- protects stainless steel instruments from irreparable damage and ensures a long-term maintenance of value
- prevents the drying of surgery residues; instruments can be stored for up to 72 hours
- facilitates subsequent cleaning
- improves personnel safety associated with reprocessing instruments

Technical data: pH-value 9.5 (undiluted, 20 °C) Density Approx.. 1.0 g/cm3 (20 °C)

Add neodisher MultiZym Instrument Cleaning Detergent w/Triple Enzyme Power, 5-Liters/Container

Product code	Туре	Size	Price 1-4	5-9	10+
C-NPS-1	PreStop Spray	750ml	\$29.50	\$28.03	\$26.55 (Bottle)
C-NPS-12 (case)	PreStop Spray	12/case	\$198.50	\$188.57	\$178.65 (Case)



Purchase online



Authorized Distributor 39 Western Hwy., West Nyack, NY 10994 www.duralinesystems.com sales@duralinesystems.com 877-561-0500